







SALUTE RDX

The **Salute RDX** mattress system provides alternating pressure therapy with micro low air loss to aid patients who are at risk of developing pressure ulcers or have existing stage 1, 2, 3, or 4 ulcers. Caregiver friendly, the Salute RDX is an affordable, therapeutic surface solution for any clinical care setting.

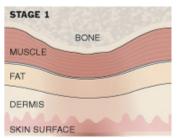

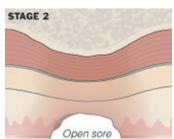

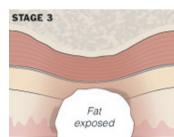

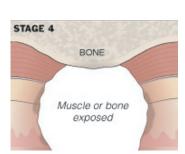



Optional Raised Bolster Cover available for added patient safety. Bolster cover has center access providing easy ingress and egress.

> Model No: Dimensions:

Weight:

MASTER CONTROL UNIT	
Dimension:	(L)13.5"x (W)7.3"x (H)8.3"
Weight:	7.1 lbs
Rating Voltage:	110-120 volts
Max Current:	0.2A max.
Rating Frequency:	60Hz
Cycle Time:	5, 10, 15, 20 mins
	Static function
	3 in 1 alternation
Min/Max Pressure:	12~47 mmHg +/- 6 mmHg
Air flow rate:	> 6 l/min
Protection Type:	Class I, Type BF
Low Pressure Alarm:	Visual / Audio
Power Failure Alarm:	Visual / Audio
	Maintenance service LED indicator

Weight Capacity:	500 lbs
Cell Material:	Black Nylon with PU coating
	Cell in cell design to stay inflated during power outage
Cell Number:	18 cells
Cover Material:	Blue Nylon with PU coating finish
Base Material:	Black Polyester with PVC coating
	CPR Quick Release System
	Cable Management System on bottom cover
HCPSC Code	E0277

Salute RDX

25.3 lbs

(L)80"x (W)36"x (H)10"

PRIUS HEALTHCARE USA

Creating innovation and difference

T (813) 854-5464
E info@priushc.com
W priushcusa.com

MATTRESS REPLACEMENT

4027 Tampa Road, #3200, OLDSMAR, FL 34677

SALUTE RDX

Provides you, the caregiver, with peace of mind, and your patient with world class care and comfort





Perfect for




Nylon covered bottom cover for friction and noise reduction when raising bed's head pan





Slimline pump design is unobstructive

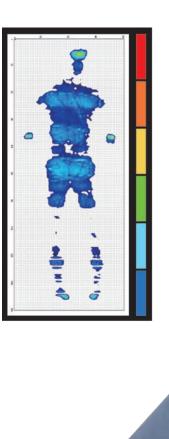

PRESSURE MAPPING

Salute RDX provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.



- Interchangeable air cell design provides easy cleaning and minimizes the proliferation of bacteria
- The Salute RDX employs 8" air bladders with cell in cell design (4.5" of therapy, 3.5" of static air), to provide therapeutic alternating pressure and also allow the mattress to stay inflated in a power outage.
- The Salute RDX replacement mattress can accommodate patient weights up to 500 lbs







3-in-1 alternating pressure covers 66% of patient's body (vs. 50% on 2 in 1 systems) to provide a more comfortable

yet still effective therapy





system deflates the

system within 20

seconds

Salute RDX is your

and treatment of

pressure ulcers

partner in prevention

footboard